




Stealth Products

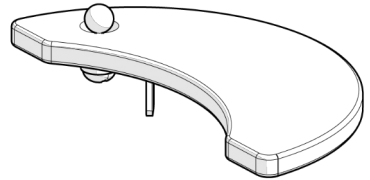

Installation Manual

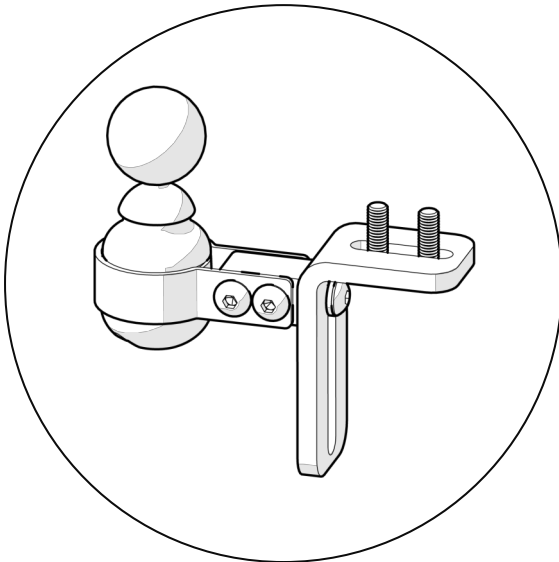

GAT440M-51M-L

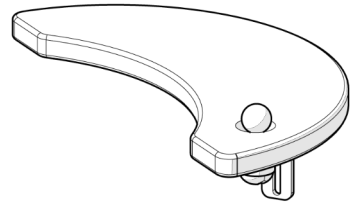


GAT440M-51M-R

GAT440M-51M (Left/Right)

Micro Joystick Assembly

Compatible GAT400 Series Mounting

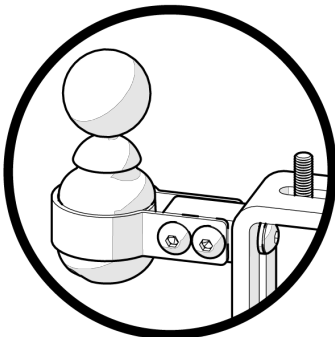

INSTALLATION INSTRUCTIONS FOR

Stealth Products: IDM-MICRO-9

mo-Vis: P002-51 Micro Joystick

Stealth Commitment

1.1

Stealth Products strives for 100% customer satisfaction. Your complete satisfaction is important. Please contact us with feedback or suggested changes that will help improve the quality and usability of our hardware. You may reach us at:



Stealth Products, LLC

104 John Kelly Drive, Burnet, TX 78611

Phone: (512) 715-9995

Fax: (512) 715-9954

Toll Free: (800) 965-9229

Toll Free: (800) 806-1225

info@stealthproducts.com

www.stealthproducts.com

EU REP PRIDE MOBILITY PRODUCTS EUROPE BV De Zwaan 3 1601 MS Enkhuizen The Netherlands	CH REP SODIMED SA Chemin Praz Devant 12 1032 Romanel sur-Lausanne Switzerland	UK Authorized Representative PRIDE MOBILITY PRODUCTS LTD. 32 Wedgwood Road, Bicester, Oxfordshire OX26 4UL United Kingdom	UK CA
---	---	---	---------------------------------------

Installation Date :

Dealer :

Dealer Stamp :

--

Serial Number :

--

Warning Labels

2.1

Warnings are included for the safety of the user, client, operator and property. Please read and understand what the signal words **SAFETY, NOTICE, CAUTION, WARNING** and **DANGER** mean, how they could affect the user, those around the user, and property.

 DANGER	Identifies an imminent situation which (if not avoided) will result in severe injury, death, and property damage.
 WARNING	Identifies a potential situation which (if not avoided) will result in severe injury, death, and property damage.
 CAUTION	Identifies a potential situation which (if not avoided) will result in minor to moderate injury, and property damage.
<i>NOTICE</i>	Identifies important information not related to injury, but possible property damage.
SAFETY	Indicates steps or instructions for safe practices, reminders of safe procedures, or important safety equipment that may be necessary.

Limited Liability

2.2

Stealth Products, LLC accepts no liability for personal injury or damage to property that may arise from the failure of the user or other persons to follow the recommendations, warnings, and instructions in this manual.

Stealth Products does not hold responsibility of final integration of final assembly of product to end user. Stealth Products is not liable for user death or injury.

Testing

2.3

Initial setup and driving should be done in an open area free of obstacles until the user is fully capable of driving safely.

The hardware should always be tested without any person sitting in the wheelchair, and until every alteration of the physical installation or adjustment of the joystick is complete.

1.0 Customer Satisfaction i

1.1 Stealth Commitment i

2.0 Warnings ii

2.1 Warning Labels ii

2.2 Limited Liability ii

2.3 Testing ii

3.0 Table of Contents..... iii

4.0 Parts and Accessories..... 1

4.1 Mounting Hardware Parts 1

4.2 Tools 1

5.0 Installation Instructions 2

5.1 Installing Mount and Joystick 2

5.2 Mounting Height Adjustment 3

6.0 Maintenance 4

6.1 Maintenance 4

6.2 Cleaning 4

6.1 Disinfecting 4


7.0 Warranty 5

7.1 In Case of Product Failure 5

Mounting Hardware Parts


4.1

GAT440M-51M includes:

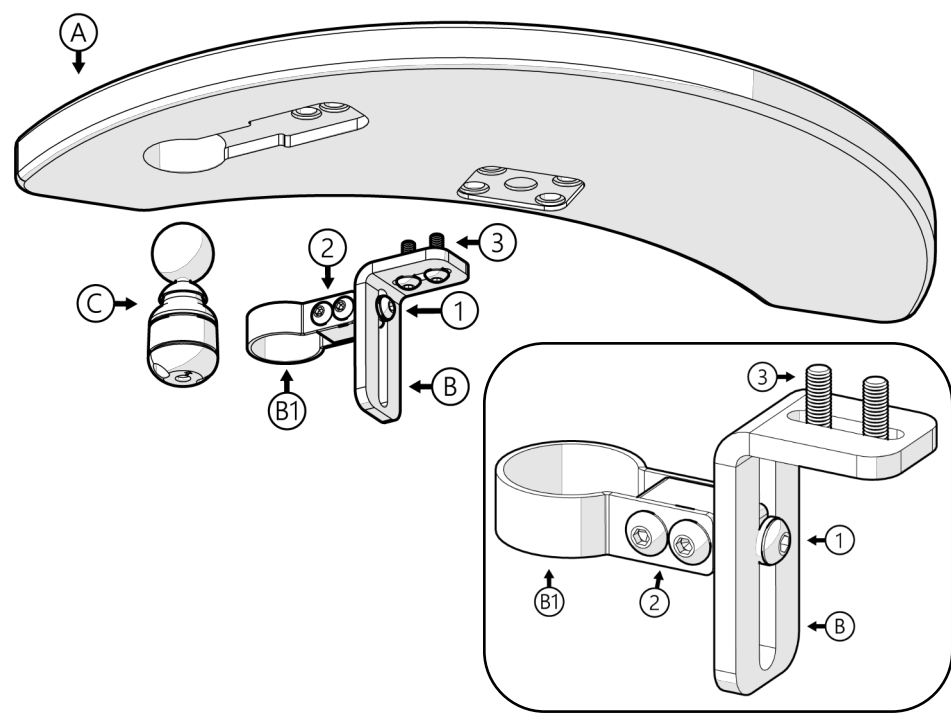

Screws& Tools

4.2

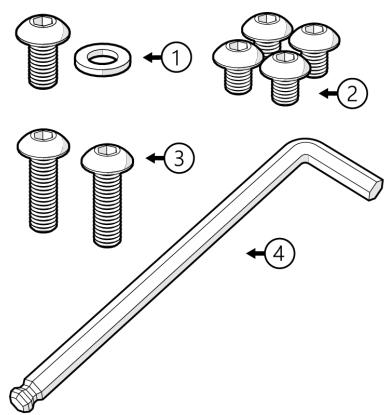


A	Eclipse/Mini Tray		
B	Mount	B1	Mounting Strap
C	Micro Joystick ¹		
1	10-32 x 1/4" BHS	1A	1/4" Steel Washer
2	10-32 x 3/8" BHS		
3	10-32 x 5/8" BHS		
4	5/32" Allen (Hex Key)		
1. Not included			

⚠ CAUTION

DO NOT OVER TIGHTEN SCREWS. This could cause wearing of screws or joystick threads, causing the client to lose efficient control.

Installing Joystick

5.1

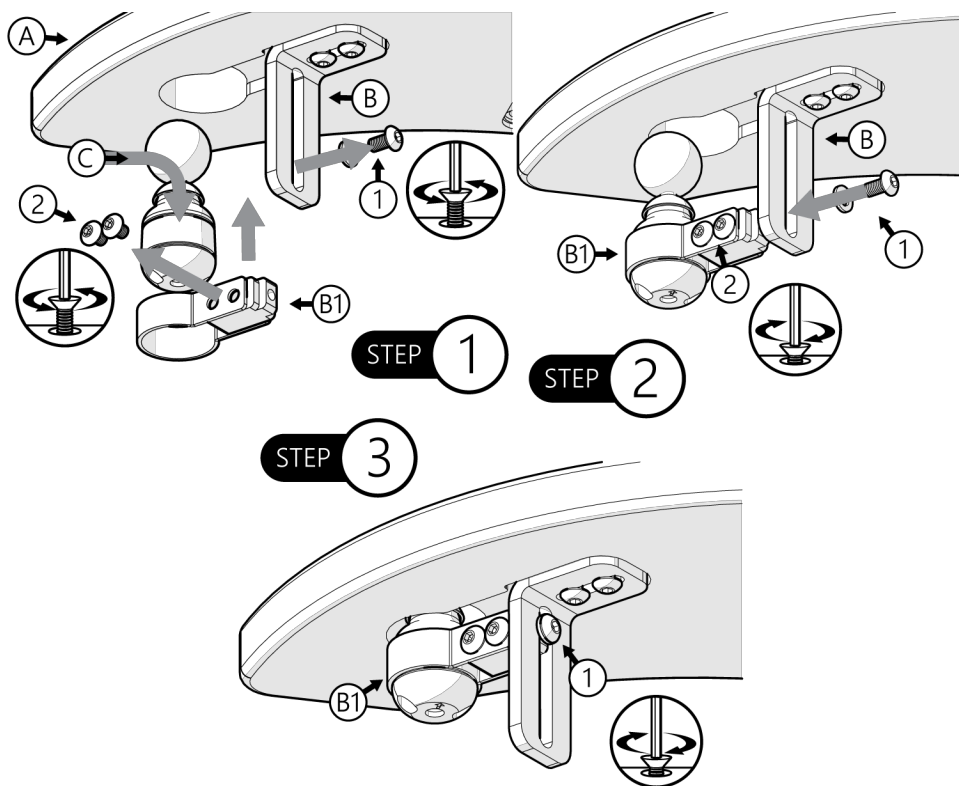


Step 1. Release mounting clamp (B1) from GAT440M-51M mount (B), loosening clamping screw (1). Loosen screws (2) from one of the sides of mounting clamp (B1) and insert Micro Joystick (C). Micro Joystick comes with default mounting clamp, remove it.

Step 2. Tighten screws (2) from mounting clamp (B1) and set into position in the mount (B) slot. Replace clamping screw.


Step 3. Tighten the screw (1) with 5/32" Allen Wrench (Hex Key) (4), but do not over tighten.

Mounting Height Adjustment

5.2

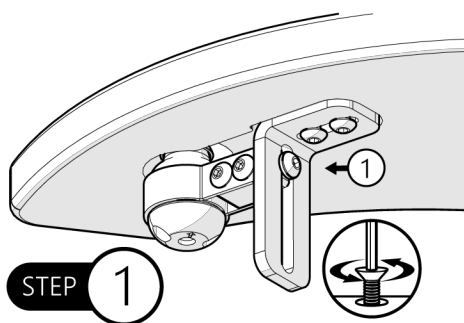

Step 1: Loosen clamping screw (1).

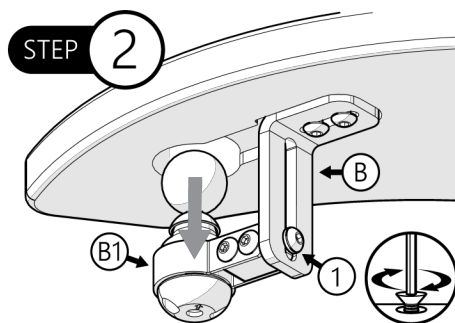

Step 2: Slide mounting clamp (B1) in mount (B) adjustment slot to the desired position. Tighten screw (1) when done.

NOTE

Mount has a 2in. (5.08 cm) range of adjustment.

SAFETY

For joystick maintenance guidelines refer to the device's user manual.

Maintenance**6.1**

- Check every 6 months that screws are tight. Gently tighten any screws that may have become loose.
- Maintain that the hardware is clean and free of dust.

Cleaning**6.2**

- Wash away surface dirt and grime with a mild dish detergent and rinse away.
- Moisten a soft cloth with water to rub away any remaining dirt.
- Rinse thoroughly with clean water.

Disinfecting**6.3**

- Use a damp cloth with any disinfectant detergent and lightly rub hardware surface.
- Rinse thoroughly with clean water.

Our products are designed, manufactured, and produced to the highest of standards. If any defect in material or workmanship is found, Stealth Products will repair or replace the product at our discretion. Any implied warranty, including the implied warranties of merchantability and fitness for a particular purpose, shall not extend beyond the duration of this warranty. Stealth Products, LLC does not warrant damage due to, but not limited to:

- Misuse, abuse, or misapplication of products.
- Modification of product without written approval from Stealth Products, LLC.
- Any alteration or lack of serial number, where applicable, will automatically void this warranty.
- Stealth Products, LLC is liable for replacement parts only.
- Stealth Products, LLC is not liable for any incurred labor costs.

No person is authorized to alter, extend, or waive the warranties of Stealth Products, LLC.

Stealth Products warrants against failure due to defective materials or workmanship:

Covers: 2 years

Hardware: 5 years

Electronics: 3 years

In Case of Product Failure

In the event of product failure covered by our warranty, please follow the procedures outlined below:

1. Call Stealth at +1 (512) 715-9995 or toll free +1-800-965-9229.
2. Request the Returns Department or obtain an RA from the Returns Department and follow department or documentation instructions.

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.



Stealth Products, LLC. • info@stealthproducts.com • www.stealthproducts.com
+1(800) 965-9229 | +1(512) 715-9995 | 104 John Kelly Drive, Burnet TX 78611