



Thoracic Support

Posture Support, H-Harness, Backpack Straps

Intended Use


The **Stealth Positioning™** line of pelvic belts, vests, chest straps and foot positioners are intended to help wheelchair users attain proper alignment of the torso, pelvis and feet. Comfort, stability and improved clinical posture are achieved through the use of breathable material and structured support.

Thoracic positioning impacts proper posture and safety. Thoracic trunk supports are placed in front of the user's torso and are designed to keep the trunk aligned with the back of the seating system. A properly positioned trunk allows for

- Improved respiratory capacity;
- Good visual field;
- Proper head positioning and neck alignment;
- Improved overall pressure distribution;
- The reduction of negative postural tendencies that can lead to orthopedic deformities and other secondary medical complications; and,
- A balance between support and function.

Installation and Use




Posture Support Shoulder Straps and Harnesses

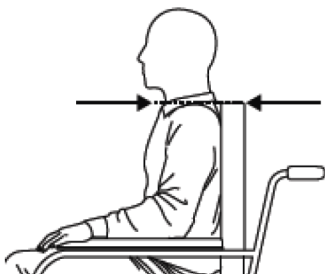

Solid seat back should be roughly at shoulder height for anterior support to fit properly. If the strap is attached lower than the shoulders, it will pull the shoulders downward. If the shoulder harness is mounted too high, it will be less effective in pulling the shoulders back.

- Shoulder straps should be positioned so the vertical straps cross mid-clavicle and the horizontal straps sit mid-trunk.
- Backpack Straps should be positioned so that vertical straps cross mid-clavicle.
- Lower straps should be installed with the included D-Rings on the back of the backrest so the strap passes under the lowest rib and does not interfere with other positioning hardware.


If the seat back height is lower than shoulder height, consider using strap risers for correct positioning. To prevent improper alignment, it is recommended that the bottom straps be adjusted first, then the top straps.

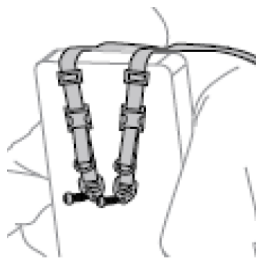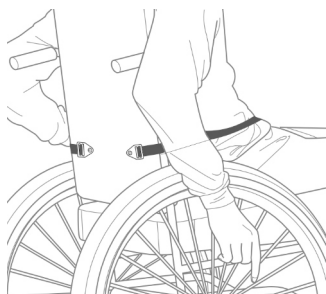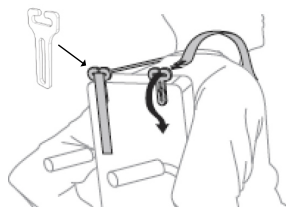

Fasten belt anchors to the seat back with the included fasteners and attach the belts to the anchors. This will allow proper adjustment of the tension straps.

Supports and harnesses are delivered in a rear-pull configuration. To switch the support or harness to a front-pull configuration, simply remove the entire strap from the harness at the tri-slide and replace it from the other end.

Install the strap section with the D-Ring and Male Buckle End to achieve a front-pull configuration. Installing the other end of the strap will result in a rear-pull configuration.

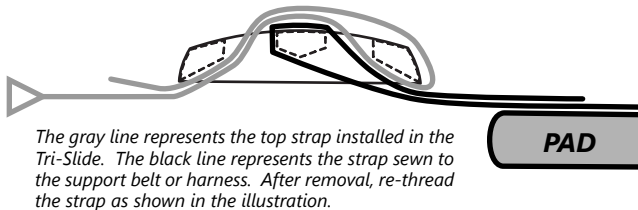

PAD

WARNING: Never use a chest harness or strap without a pelvic belt. There is a risk of chest compression or choking if the user slips down in the chair.

Swivel Buckle

The Swivel Buckle is used on these **Stealth Positioning™** products:


Positioning Systems

Y-Straps

Step 1: To release the swivel buckle from the belt, push down on the button in the center of the male swivel buckle.

Step 2: Remove the male swivel buckle and attached belt from the female swivel buckle. Make any adjustments necessary.

Step 3: To reattach the swivel buckle, slide the buckle into the housing. An audible click should occur when the male buckle has been fastened into place, indicating that the belt is secure. *If the button does not latch and hold secure, immediately refrain from using the product.*

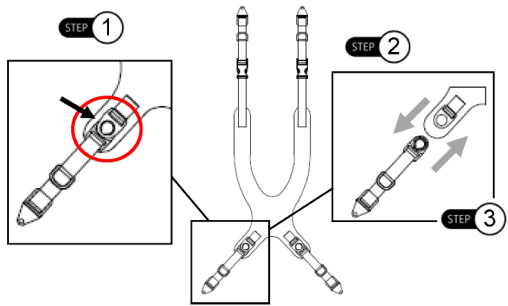

Side-Release Buckle

The Side-Release Buckle is used on these **Stealth Positioning™** products:

Pelvic Belts
Chest Straps

Backpack Straps
Foot Positioners


Y-Straps
Positioning Systems

Shoe Holders

Step 1: To remove the side-release buckle from the belt, pinch the buttons on the sides of the buckle. This will disengage the locking mechanism.

Step 2: Slide the male portion of the buckle from the housing. Make any adjustments necessary.

Step 3: To reattach the side-release buckle, slide the male portion of the buckle into the housing. An audible 'click' should occur, indicating that the belt is secure. *If the buckle does not latch and hold, immediately refrain from using the product.*

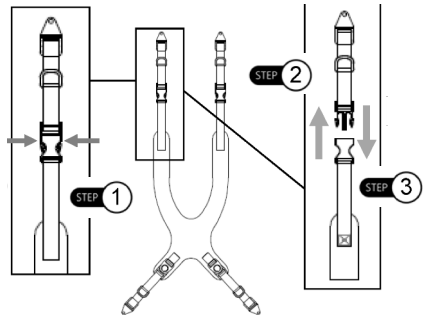

Care and Cleaning

We recommend that all products be examined on a regular basis for wear and function. Check the stitching on the webbing for stress or damage. If loose, frayed or broken parts are discovered, discontinue use and contact your provider. Repairs should only be performed by a wheelchair service professional and with genuine **Stealth Products** parts, fasteners or items appropriate for the service required.

Surfaces may be cleaned with a damp cloth and a mild, non-abrasive cleaner. Ensure all cleaners are approved for finished steel, aluminum, plastic and upholstered surfaces. Vests and straps that are removable and DO NOT contain metal components may be machine washed warm (40°C), permanent press cycle, mild process.

DO NOT bleach. DO NOT tumble dry. Dry flat. DO NOT iron.

Important

Never use a chest harness without a positioning belt. Positioning belts and supports should only be used to help support the user's posture. The belts, posture supports and straps must be snug, but must not interfere with the user's breathing. Check to make sure the user cannot slide down in the seating system, as they may suffer chest compression and may suffocate due to pressure from the belts or harnesses. A seat with angle or built-in contours may help prevent the user from sliding in the seat. Due to the possibility of a user sliding while in-chair, unsupervised use of anterior thoracic support devices could pose a serious choking hazard. Continually ensure that the user is sitting upright and that the thoracic support device does not obstruct the user's breathing in any way.

These devices do not conform to ANSI/RESNA WC20 or SAE J2249 and are not designed or approved to provide restraint or protection in the event of a crash. Wheelchair users should transfer to a vehicle seat and use the vehicle seat belt system or a child safety seat that complies with federal safety standards. The wheelchair should then be stored and secured in the vehicle. Children weighing less than 48 lbs. must always transfer to an approved child safety seat.

Before you install or begin using any product from **Stealth Products' Stealth Positioning™** line, it is important that you read and understand the content of these installation and operating instructions. Product misuse due to failure to adhere to instructions and warnings in this document will void your warranty and may result in property damage, injury or death. Always keep the operating instructions in a safe place so they may be referenced as necessary. If this document contains information you do not understand or if there are concerns about safety or operation, contact your provider.

All information, pictures, illustrations and specifications are based on the product information available at the time of printing, and are representative examples. They are not intended to be exact descriptions of the various parts of the product. Instructions are written with the expressed intent of use with standard configurations. They also contain important safety and maintenance information. For assistance with advanced configurations or options, contact your provider or **Stealth Products** at (512) 715-9995 or toll free at (800) 965-9229.

Immediately discontinue use if any function is compromised, if parts are missing or loose or if any component shows signs of excessive wear. Consult with your provider for repair, adjustment or replacement. You can download additional copies of this manual by visiting <https://stlpro.site/stealth-docs/> and searching for "Stealth Positioning" in the document manager.

Warranty

Our products are designed, manufactured, and produced to the highest of standards. If any defect in material or workmanship is found, **Stealth Products** will repair or replace the product at our discretion. Any implied warranty, including the implied warranties of merchantability and fitness for a particular purpose, shall not extend beyond the duration of this warranty. **Stealth Products, LLC** does not warrant damage due to, but not limited to, misuse, abuse, or misapplication of products and modification of product without written approval from **Stealth Products, LLC**.

Any alteration or lack of serial number, where applicable, will automatically void this warranty. **Stealth Products, LLC** is liable for replacement parts only. **Stealth Products, LLC** is not liable for any incurred labor costs. No person is authorized to alter, extend, or waive the warranties of **Stealth Products, LLC**.

To learn more about our warranty and returns policy, please visit <https://stlpro.site/returns-n-warranty>.

About Us

Stealth Products is committed to 100% customer satisfaction. Your complete satisfaction is important to us. Please contact us with suggestions or feedback to help us improve the quality and usability of our products.

You may reach us at:



Stealth Products, LLC

104 John Kelly Drive,
Burnet, TX 78611

Phone: (512) 715-9995 Toll Free: (800) 965-9229

Fax: (512) 715-9954 Toll Free: (800) 806-1225

info@stealthproducts.com - www.stealthproducts.com

EU REP
PRIDE MOBILITY PRODUCTS EUROPE BV.
De Zwaan 3, 1601 MS Enhuizen, The Netherlands

CH REP
SODIMED SA
Chemin Praz Devant 12 1032 Romanel sur-Lausanne, Switzerland

UK Authorized Representative
PRIDE MOBILITY PRODUCTS LTD.
32 Wedgwood Road Bicester, Oxfordshire OX26 4UL, United Kingdom

