

### TECHNICAL DATA

#### INCLUDES:

Comfort Plus Headrest Pad  
1" Mounting Ring

#### CB325

Combo Series Headrest Pad





#### HMO400

Headrest Mounting Hardware

### DIMENSIONS & MOUNTING

#### MOUNTING:

The Combo Series Headrest attaches to a 17mm ball mount via a SD-Ring. When you have the mounting ring in place, screw it into the adjacent holes on the headrest.

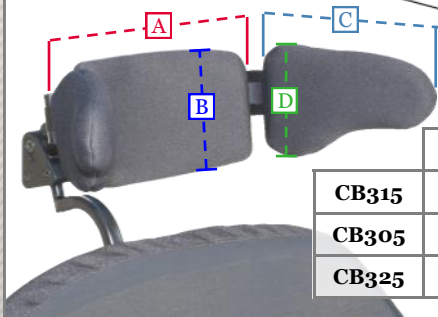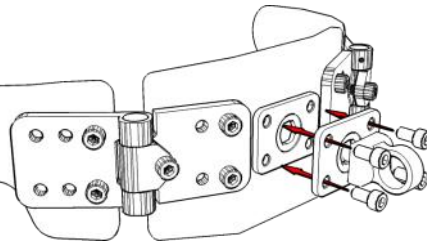

	Size	Width	Height	Width	Height
<b>CB315</b>	Small	5"	3"	3 <sup>1</sup> / <sub>8</sub> "	3 <sup>3</sup> / <sub>4</sub> "
<b>CB305</b>	Large	5 <sup>3</sup> / <sub>4</sub> "	3 <sup>1</sup> / <sub>2</sub> "	4 <sup>3</sup> / <sub>4</sub> "	3 <sup>1</sup> / <sub>2</sub> "
<b>CB325</b>	Extra Large	5 <sup>3</sup> / <sub>4</sub> "	6 <sup>1</sup> / <sub>4</sub> "	7"	6 <sup>1</sup> / <sub>4</sub> "

### PRODUCT DESCRIPTION

The Combo Series is an adjustable-curved head positioning system. It's intuitive function lets you set the shape, then reliably holds your settings. The Combo Series head positioning system consist of two lateral "wing" pads. The lateral pads were designed to provide a high degree of pivoting.

### TECHNICAL DATA

1. 5/32 tool for socket-head and flat-head screws on mounting ring and mounting hardware.
2. Includes: ball mounting ring and headrest pad.
3. Headrest attaches to: HMO400, HMO475, TWB480-1, TWB480-FDM, and TWB485-1.

### SELLING POINTS

1. The high degree of pivoting movement accommodates a wide variety of head positioning needs.
2. The built in facial laterals of the Combo Series provides amazing positioning.
3. Compatibility with all of our standard mounting hardware.
4. Ball mount allows for angle and rotational adjustment.

### BUILD YOUR COMFORT PLUS

#### STEP 1: PAD SIZE

The Combo Series comes with several types of foam and covers. According to the individual's needs, determine the necessary pad size, whether foam or gel would be beneficial, and cover material. (See Headrest Pads on back)

CB305 — CB315 — CB325

#### STEP 2: SWING-AWAYS (SUS9 SERIES)

Facial lateral Swing-Aways for pads, switches, speaker sites and chin control devices can be added to the Combo Series. (See Swing-Away Options on back)

SUS9-B-HW — SUS9DL-DB — SUNS9-B-HW

#### STEP 3: HARDWARE

The final step is deciding what type of hardware to use for mounting the Combo Series to your chair. (See Hardware Options on back)

#### FOAM/GEL

Open Cell Foam  
Closed Cell Foam  
Action Gel

#### COVERS

Dartex  
Reverse Dartex  
Vinyl  
Sheep Skin  
Neoprene

#### HARDWARE

HMO400  
HMO475  
TWB480  
TWB480-FDM  
TWB485

### HARDWARE OPTIONS



#### **HMO400**

The HMO400 is our straight mounting hardware can be adjusted fore and aft.



#### **HMO475**

The HMO475 allows fore and aft adjustment as well as rotational adjustment



#### **TWB480-1**

The TWB480 has a rotating triangle and links that can be adjusted to any position.



#### **TWB485-1**

The TWB485 combines the links of the TWB480 and curved bar of the HMO475.



#### **TWB480-FDM**

An easy transfer flip down mount and links make Stealth's TWB480-FDM.



#### **TD101**

The Tone Deflector is design to absorb quick impact and avert slanted angle loading.

### HEADREST PADS



#### **CB315**

**Small**

Combo Series Headrest



#### **CB305**

**Large**

Combo Series Headrest



#### **CB325**

**Extra Large**

Combo Series Headrest

### SWING-AWAY OPTIONS

More information available at [www.stealthproducts.com](http://www.stealthproducts.com)



#### **SUS9-B-HW**

Stealth Product's Swing-Away hardware offers a sleek, unobtrusive design while featuring an easy, one-hand release mechanism for swinging pads and switches back a full 100°.



#### **SUS9DL-DB**

The dual swing away design allows for more attachment sites for facial pads or switches for those clients who need more real estate or extra support.



#### **SUNS9-B-HW**

Convenient sizing, versatility, comfort and intuitive function with a clean, elegant design makes Stealth's SUNS9 the best choice when in need for pediatric facial pads and switches.



#### **900**

**Facial Pad**  
(4.5" x 2")



#### **910**

**Facial Spot Pad**  
(3" x 2")



#### **920**

**Long Facial pad**  
(6" x 2")



#### **930**

**Long Niño Facial Pad**  
(6" x 1.5")



#### **F8675**

Stealth's figure 8 bracket adds real estate for an extra attachment.



#### **JOYSTICK MOUNTS**

**SM600:** Mini Proportional mount  
**SM610:** 1812 mount  
**SM620:** Tash Mini mount  
**SM630:** ASL MEC mount



#### **SP625**

Sip & Puff mount



#### **CN600**

Swing-away flex Micro-Lite switch.