



**STEALTH
POSITIONING**

stealthproducts.com

***stability in a
Snap!***

Stealth Products
Solutions For Life®





The Stealth Positioning™ line of pelvic, trunk, and lower extremity supports is designed to enhance stability, support proper alignment, and encourage healthy posture for wheelchair users. Engineered for quick application and a secure fit, these elements help increase the user's ability to participate in daily activities. With Stealth Products®, stability is a snap!



Breathable Material

Constructed from a unique open-cell material that allows air to circulate and moisture to escape. This contributes to reduced heat buildup and a cooler, drier, and more comfortable user experience.



Offset Belt Anchors

A slight offset allows for easy adjustment without loosening the anchor bolts.



Swivel Buckles

Swivel buckles on the lower attaching straps allow self-adjustment, supporting the correct angle of attachment.



Dynamic Support

Dynamic Support provides postural control, a degree of trunk movement, and a softer fit on the body. The user has movement within a limited range encouraging return to their neutral position.



Structured Support

The structured material is more resistant to high tone which helps maintain correct posture.



Full Length Attaching Straps

Provide for proper positioning of the top side release buckles, ensuring that the posture support maintains adjustment.

Posture Support

- Breathable Material
- Dynamic & Structured Support
- Classic & Contoured Styles
- Swivel Buckles
- Full Range of Sizes
- Offset Belt Anchor
- Large Neck Cavity
- Full Length Attaching Straps

CLASSIC

The Classic Style offers a large padded surface for the thoracic area.




	A	B	C	D	E	F	G	H
XS	7.5	12.5	2	3.5	7.5	3.5	3	2
S	8.5	15	2	4.5	9	4	4	2
M	9.5	19	2	5.5	11	5	5	2.75
L	10.5	22.5	2	6.5	13	6	6	3.5
XL	11.5	25.5	2.5	6.5	14	7	7	4.5

	XS	S	M	L	XL
Dynamic	IFPS1101	IFPS1102	IFPS1103	IFPS1104	IFPS1105
Structure	IFPS1201	IFPS1202	IFPS1203	IFPS1204	IFPS1205

CONTOURED

The Contoured Style has a narrowed section over the sternum that allows a better fit for females, G-tube users or those requiring a smaller support surface.

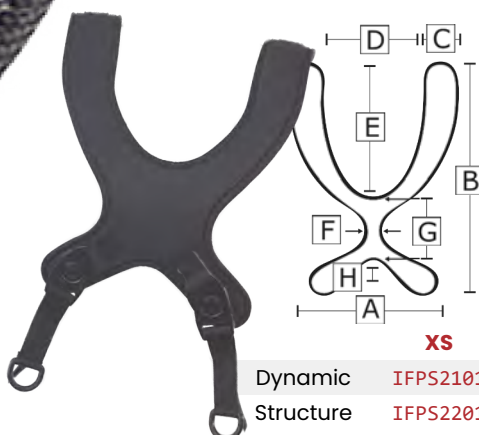

	A	B	C	D	E	F	G	H
XS	7.5	12.5	2	3.5	7.5	1	3.5	1.25
S	8.5	15	2	4.5	9	1	4.25	1.75
M	9.5	18.75	2	5.5	11	1	5.25	2.5
L	10.5	22.5	2	6.5	13	1.25	6	3.5
XL	11.5	25.5	2.5	6.5	14	1.5	7	4.5

	XS	S	M	L	XL
Dynamic	IFPS2101	IFPS2102	IFPS2103	IFPS2104	IFPS2105
Structure	IFPS2201	IFPS2202	IFPS2203	IFPS2204	IFPS2205



CLASSIC VINYL FACE (DYNAMIC)

IFPS1401	Pull Down, Extra Small (XS)
IFPS1402	Pull Down, Small (S)
IFPS1403	Pull Down, Medium (M)
IFPS1404	Pull Down, Large (L)
IFPS1405	Pull Down, Extra Large (XL)

CLASSIC VINYL FACE (STRUCTURED)

IFPS1301	Pull Down, Extra Small (XS)
IFPS1302	Pull Down, Small (S)
IFPS1303	Pull Down, Medium (M) •
IFPS1304	Pull Down, Large (L)
IFPS1305	Pull Down, Extra Large (XL)



CONTOURED VINYL FACE (DYNAMIC)

IFPS2401	Pull Down, Extra Small (XS)
IFPS2402	Pull Down, Small (S)
IFPS2403	Pull Down, Medium (M) •
IFPS2404	Pull Down, Large (L)
IFPS2405	Pull Down, Extra Large (XL)

CONTOURED VINYL FACE (STRUCTURED)

IFPS2301	Pull Down, Extra Small (XS)
IFPS2302	Pull Down, Small (S)
IFPS2303	Pull Down, Medium (M)
IFPS2304	Pull Down, Large (L)
IFPS2305	Pull Down, Extra Large (XL)



CONTOURED ZIPPER FACE (DYNAMIC)

IFPS3102	Pull Down, Small (S)
IFPS3103	Pull Down, Medium (M)
IFPS3104	Pull Down, Large (L)
IFPS3105	Pull Down, Extra Large (XL)

CLASSIC ZIPPER FACE (STRUCTURED)

IFPS3202	Pull Down, Small (S)
IFPS3203	Pull Down, Medium (M)
IFPS3204	Pull Down, Large (L) •
IFPS3205	Pull Down, Extra Large (XL)

BACKPACK, PAIR (STRUCTURED)

IFBP4201	Pull Down, Extra Small (XS)
IFBP4202	Pull Down, Small (S)
IFBP4203	Pull Down, Medium (M) •
IFBP4204	Pull Down, Large (L)






H HARNESS (STRUCTURED)

IFHH5201	Pull Down, Extra Small (XS)
IFHH4202	Pull Down, Small (S)
IFHH4203	Pull Down, Medium (M) •
IFHH4204	Pull Down, Large (L)



CHEST STRAPS ALL CHEST STRAPS ARE DYNAMIC

Y STRAPS

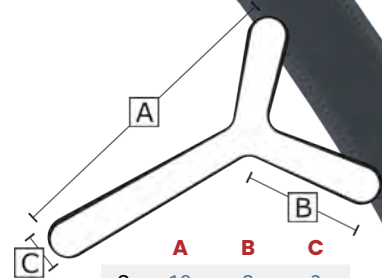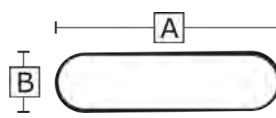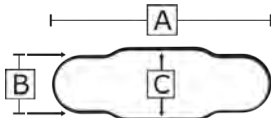

		A	B	C
Belt Anchor	XS	14	2.8	3.5
	S	16	2.8	3.75
	M	18	3.1	4
	L	21	3.1	4.5
	XL	23	3.25	4.75

		A	B
Hook & Loop	XS	7	2
	S	8.5	3
	M	10.25	4
	L	11.5	5
	XL	13	6

	XS	S	M	L	XL	Non-Padded
Belt Anchor	IFCS4111	IFCS4112	IFCS4113	IFCS4114	IFCS4115	IFCS4221
Hook & Loop	IFCS4121	IFCS4122	IFCS4123	IFCS4124	IFCS4125	IFCS4222


	S	M	L
Dynamic	IFYS3102	IFYS3103	IFYS3104
Structure	IFYS3202	IFYS3203	IFYS3204



Pelvic Belts & Pads


- Help maintain a neutral pelvic position, which directly affects seated body alignment and distributes pressure appropriately.
- Increase postural stability, decreasing the effort needed to sit upright and facilitating upper extremity movement.
- Contribute to proper head positioning and neck alignment.
- Reduce negative postural tendencies that can lead to orthopedic asymmetries and other secondary medical complications.

2 POINT PAD SIZES

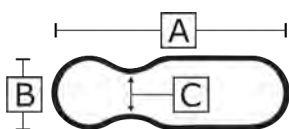


	Part Number	A	B	C
XS	IFPB101	4.25	1.75	1.25
S	IFPB102	5.5	2.25	1.25
M	IFPB103	7.5	2.5	1.25
L	IFPB104	9.5	2.25	2.25
XL	IFPB105	11.5	2.78	2.25

4 POINT PAD SIZES

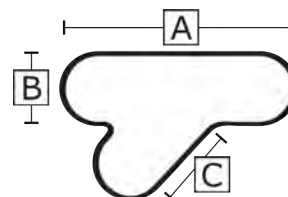


	A	B	C
XS	4.25	1.75	3.25
S	5.5	2.25	4
M	7.5	2.5	4
L	9.5	2.25	4

BELT SIZES

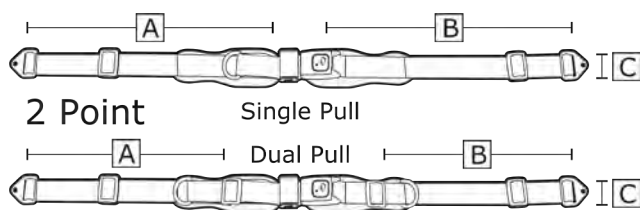

Style	A	B	C	Buckle	Pad Sizes
Single Pull	32	27	1	Push Button	NP / XS
Single Pull	32	27	1	Side Release	NP
Dual Pull	24	24	1	Push Button	XS / S
Single Pull	32	27	1.5	PB / SR	NP / M / L
Dual Pull	24	24	1.5	PB / SR	M / L
Single Pull	32	27	2	Push Button	NP / M / L
Single Pull	32	27	2	Side Release	L / XL
Dual Pull	24	24	2	PB / SR	L / XL

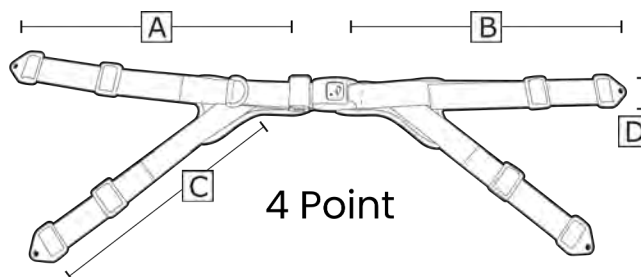

A	B	C	D	Buckle	Sizes
32	27	20	1	Push Button	XS
32	27	20	1.5	Push Button	S / M / L

Leg Harness



Sizes	A	B	C	D	E
Small	1	9	1.5	1	5
Medium	1	12	1.5	1.5	5
Large	1.5	15	2		

PELVIC BELTS



2 POINT NON-PADDED BELTS

PUSH BUTTON

IFPB5110	Single Pull, 1
IFPB5210	Single Pull, 1.5
IFPB5310	Single Pull, 2

SIDE RELEASE

IFPB6110	Single Pull, 1
IFPB6210	Single Pull, 1.5 •



2 POINT PADDED BELTS (DUAL PULL)

1" BELTS

IFPB5121	Push Button, Extra Small Pads (4.5)
IFPB5122	Push Button, Small Pads (5.5)

1.5" BELTS

IFPB5223	Push Button, Medium Pads (7.5)
IFPB5224	Push Button, Large Pads (9.5)
IFPB6223	Side Release, Medium Pads (7.5)
IFPB6224	Side Release, Large Pads (9.5)

2" BELTS

IFPB5324	Push Button, Large Pads (9.5)
IFPB5325	Push Button, Extra Large Pads (11.5)
IFPB6324	Side Release, Large Pads (9.5) •
IFPB6325	Side Release, Extra Large Pads (11.5)



4 POINT BELTS (V-CONTOUR STYLE)

1" BELTS

IFPB5141	Single Pull, Push Button, Extra Small Pads (4.5)
IFPB5242	Single Pull, Push Button, Small Pads (5.5)
IFPB5151	Rear Pull, Push Button, Extra Small Pads (4.5)
IFPB5152	Rear Pull, Push Button, Small Pads (5.5)

1.5" BELTS

IFPB5243	Single Pull, Push Button, Medium Pads (7.5)
IFPB5244	Single Pull, Push Button, Large Pads (9.5) •
IFPB5253	Rear Pull, Push Button, Medium Pads (7.5)
IFPB5254	Rear Pull, Push Button, Large Pads (9.5)
IFPB6253	Rear Pull, Side Release, Medium Pads (7.5)
IFPB6254	Rear Pull, Side Release, Large Pads (9.5)

2" BELTS

IFPB5253	Rear Pull, Push Button, Medium Pads (7.5)
IFPB5254	Rear Pull, Push Button, Large Pads (9.5)
IFPB5253	Rear Pull, Side Release, Medium Pads (7.5)
IFPB5254	Rear Pull, Side Release, Large Pads (9.5)



RETRACTABLE PELVIC BELT

IFRB5350	Pelvic Belt, Push Button, Retractable 2" Webbing, Non-Padded
----------	--



2 POINT PADDED BELTS

1" BELTS

IFPB5111	Single Pull, Push Button, Extra Small Pads (4.5)
IFPB5112	Single Pull, Push Button, Small Pads (5.5)
IFPB5131	Rear Pull, Push Button, Extra Small Pads (4.5)
IFPB5132	Rear Pull, Push Button, Small Pads (5.5)

1.5" BELTS

IFPB5213	Single Pull, Push Button, Medium Pads (7.5)
IFPB5214	Single Pull, Push Button, Large Pads (9.5)
IFPB5233	Rear Pull, Push Button, Medium Pads (7.5)
IFPB5234	Rear Pull, Push Button, Large Pads (9.5)
IFPB6213	Single Pull, Side Release, Medium Pads (7.5)
IFPB6214	Single Pull, Side Release, Large Pads (9.5) •
IFPB6233	Rear Pull, Side Release, Medium Pads (7.5)
IFPB6234	Rear Pull, Side Release, Large Pads (9.5)

2" BELTS

IFPB5314	Single Pull, Push Button, Large Pads (9.5)
IFPB5315	Single Pull, Push Button, Extra Large Pads (11.5)
IFPB5334	Rear Pull, Push Button, Large Pads (9.5)
IFPB5335	Rear Pull, Push Button, Extra Large Pads (11.5)
IFPB6314	Single Pull, Side Release, Large Pads (9.5)
IFPB6315	Single Pull, Side Release, Extra Large Pads (11.5)
IFPB6334	Rear Pull, Push Button, Large Pads (9.5)
IFPB6335	Rear Pull, Push Button, Extra Large Pads (11.5)



4 POINT BELTS (LINEAR STYLE)

1" BELTS

IFPB5145	Push Button, Extra Small Pads (4.5)
----------	-------------------------------------

1.5" BELTS

IFPB5146	Push Button, Small Pads (5.5)
IFPB5147	Push Button, Medium Pads (7.5)
IFPB5148	Push Button, Large Pads (9.5) •


LEG HARNESS



IFLH1111-S	Dynamic, 1.5 Webbing, Small (S)
IFLH1112-M	Dynamic, 1.5 Webbing, Medium (M)
IFLH1123-L	Dynamic, 2 Webbing, Large (L) •

Foot Positioning

ANKLE SUPPORT SIZES

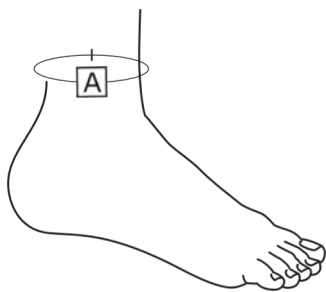

	Ankle Strap A (Circumference)	Foot Plates B (Width)	C (Length)
XS	6-8		
S	7-8.5	4.75	7.75
M	8-11	5.25	9.25
L	9-12	5.75	10.25
XL	12.5-15		



FIXED ANKLE SUPPORTS

IFFP3210- [L/R]	Structured ¾ in Webbing, Extra Small (XS)
IFFP3211- [L/R]	Structured ¾ in Webbing, Small (S)
IFFP3212- [L/R]	Structured ¾ in Webbing, Medium (M) •
IFFP3213- [L/R]	Structured ¾ in Webbing, Large (L)
IFFP3214- [L/R]	Structured ¾ in Webbing, Extra Large (XL)

FIXED ANKLE SUPPORTS (HOOK & LOOP)

IFFP3310- [L/R]	Structured H&L, ¾ in Webbing, Extra Small (XS)
IFFP3311- [L/R]	Structured H&L, ¾ in Webbing, Small (S)
IFFP3312- [L/R]	Structured H&L, ¾ in Webbing, Medium (M)
IFFP3313- [L/R]	Structured H&L, ¾ in Webbing, Large (L)
IFFP3314- [L/R]	Structured H&L, ¾ in Webbing, Extra Large (XL)



ANKLE POSITIONING SYSTEMS

IFFP1110- [L/R]	Dynamic ¾ in Webbing, Extra Small (XS)
IFFP1111- [L/R]	Dynamic ¾ in Webbing, Small (S)
IFFP1122- [L/R]	Dynamic 1 in Webbing, Medium (M) •
IFFP1123- [L/R]	Dynamic 1 in Webbing, Large (L)
IFFP1124- [L/R]	Dynamic 1 in Webbing, Extra Large (XL)



FULL FOOT POSITIONING SYSTEMS

IFFP2110- [L/R]	Dynamic ¾ in Webbing, Extra Small (XS)
IFFP2111- [L/R]	Dynamic ¾ in Webbing, Small (S)
IFFP2122- [L/R]	Dynamic 1 in Webbing, Medium (M) •
IFFP2123- [L/R]	Dynamic 1 in Webbing, Large (L)
IFFP2124- [L/R]	Dynamic 1 in Webbing, Extra Large (XL)

SHOE HOLDER SIZES



	A	B
S	5	3
M	6.25	3.5
L	7.75	4.25

SHOE HOLDERS

IFSH4211	Padded Straps ¾ in Webbing, Small (S) •
IFSH4212	Padded Straps ¾ in Webbing, Medium (M)
IFSH4213	Padded Straps ¾ in Webbing, Large (L)
IFTS4211	Padded Toe Strap Only, Small (S)
IFTS4212	Padded Toe Strap Only, Medium (M)
IFTS4213	Padded Toe Strap Only, Large (L)

PADDED FOOT SUPPORT



Moldable Upper Support
See Owner's Manual for details



* Dual always .5
narrower than
described to fit
between the leg
rest frame.



		Width A	Depth B	Standard	Tall	For Dual with divider only	
				Height C		Divider Width D	Leg Trough E
Dual* (divided or not divided)	SINGLE	7.5	10	8.5	12.5	-	-
	12	11.5	10.5	8.5	-	-	-
	14	13.5	10.5	8.5	-	-	-
	16	15.5	10.5	8.5	12.5	2.5	5
	18	17.5	10.5	8.5	12.5	2.5	6
	20	19.5	10.5	8.5	12.5	2.5	7
	22	21.5	10.5	8.5	12.5	2.5	8
	24	23.5	10.5	8.5	12.5	2.5	9

SINGLE PADDED FOOT SUPPORT

IFPFB-SS-1	Standard •
IFPFB-ST-1	Tall •

DUAL PADDED FOOT SUPPORT - TALL -

NO DIVIDER	
IFPFB-FT-16	No Divider Tall 16 Wide
IFPFB-FT-18	No Divider Tall 18 Wide
IFPFB-FT-20	No Divider Tall 20 Wide
IFPFB-FT-22	No Divider Tall 22 Wide
IFPFB-FT-24	No Divider Tall 24 Wide
WITH DIVIDER	
IFPFB-DT-16	With Divider Tall 16 Wide
IFPFB-DT-18	With Divider Tall 18 Wide
IFPFB-DT-20	With Divider Tall 20 Wide
IFPFB-DT-22	With Divider Tall 22 Wide
IFPFB-DT-24	With Divider Tall 24 Wide

DUAL PADDED FOOT SUPPORT - STANDARD -

NO DIVIDER	
IFPFB-FS-12	No Divider Standard 12 Wide
IFPFB-FS-14	No Divider Standard 14 Wide
IFPFB-FS-16	No Divider Standard 16 Wide
IFPFB-FS-18	No Divider Standard 18 Wide
IFPFB-FS-20	No Divider Standard 20 Wide •
IFPFB-FS-22	No Divider Standard 22 Wide
IFPFB-FS-24	No Divider Standard 24 Wide
WITH DIVIDER	
IFPFB-DS-16	With Divider Standard 16 Wide
IFPFB-DS-18	With Divider Standard 18 Wide
IFPFB-DS-20	With Divider Standard 20 Wide •
IFPFB-DS-22	With Divider Standard 22 Wide
IFPFB-DS-24	With Divider Standard 24 Wide

Calf Support

CALF STRAP




CALF STRAP PART NUMBERS

(STRUCTURE)

IFCF1511-S	Small (S)				
IFCF1512-M	Medium (M)				
IFCF1513-L	Large (L)				
	Sizes	A	B	C	D
	Small	31.5	14	4.375	3.5
	Medium	35.5	16	4.875	3.5
	Large	39.5	18	5.375	3.5

CALF PANEL

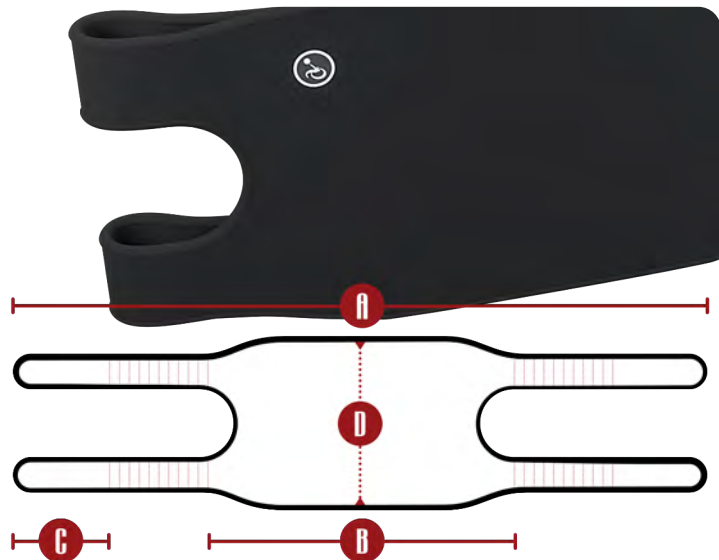

CALF PANEL PART NUMBERS

(STRUCTURE)

IFCF2511-S	Small (S)			
IFCF2512-M	Medium (M)			
IFCF2513-L	Large (L)			
Sizes	A	B	C	D
Small	31.5	14	4.375	7.5
Medium	35.5	16	4.875	8.25
Large	39.5	18	5.375	9.25



Accessories



BUCKLE COVER

IFMDBC10 1" Buckle
IFMDBC15 1.5" Buckle



CAM LOCK

IFMDCM10 1" Webbing
IFMDCM15 1.5" Webbing
IFMDCM20 2" Webbing



CANE CLAMP KIT

IFMDCC Pair



Stealth Products
Solutions For Life®

104 John Kelly Drive, Building A
Burnet, Texas 78611

(800) 965-9229

stealthproducts.com